


# Auto-disable Syringe



## Features

- Pull back the plunger to charge the syringe with solution.
- Press plunger forward to complete injection until it reaches the stop position. Lock mechanism will be activated and locks plunger in stop position. An audible "Click" will sound.
- Forcing the plunger backward will cause it to break dispose safety in a disposable container.

## Operation process

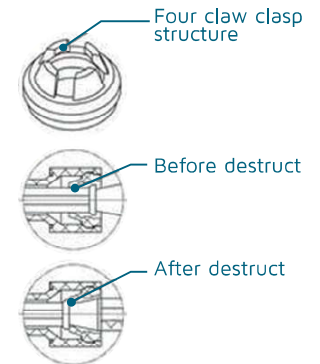

Remove protect cap



Pumping drugs



Pull plunger until stop ring



Complete self-destruct after drugs injection

## Specifications

Syringe Capacity	Needle	
	Gauge	Length
0.1ml	25G	5/8"(16mm)
0.5ml	25G	3/4"(19mm)
1.0ml	25G	3/4"(19mm)
2.0ml	25G	1.0"(25mm)
2.5ml	25G	1.0"(25mm)
3.0ml	25G	1.0"(25mm)
5.0ml	25G	1.0"(25mm)

▲ Oem & odm are available on request.